




MoleScope™ II

Quick Start Guide for A/A+/B/B+ Models

P002 - 004 - QSGB
2023-06-06 V4.2

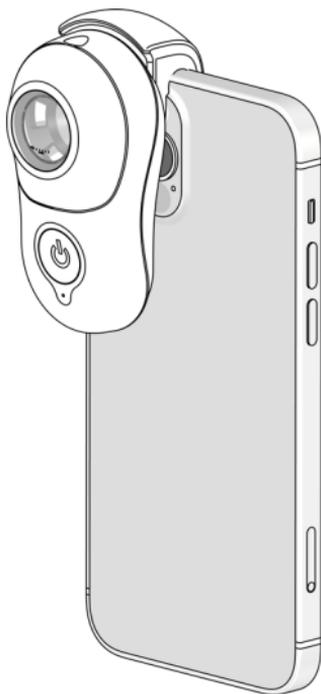

Consult Instruction for Use
www.molescope.com/molescope2




Follow instructions for use



Step 1:

Install the rubber insert as shown in the image. Make sure the rubber insert is pushed in and is tight in place.

Your device may come with a pre-installed rubber insert. If so, please skip step 1.

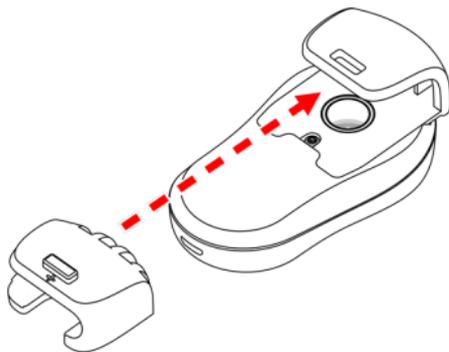


Rubber Insert

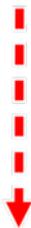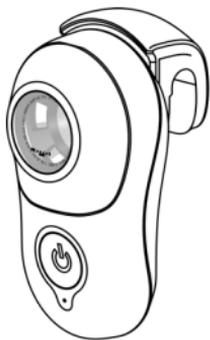

Step 2:

Slide the device over your phone camera and attach tightly. You may slide it from top or side based on you phone model. Make sure your device lens is centered over your smartphone camera.

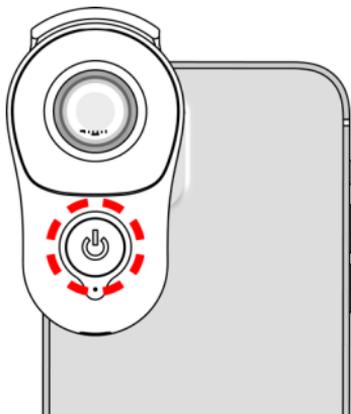

Step 3:

Turn on the light prior your imaging session by pressing the power button.

Step 4:

Gently press the MoleScope flat against the skin making sure the lesion is centered. Tap your mobile device screen to focus and take the image.



Step 5:

Your MoleScope comes with a contact cap and there is an extra cap in the box for non-contact imaging. Simply lift up the the contact cap and replace with the non-contact cap if needed.

Make the most of MoleScope II with DermEngine


For medical professionals and care providers, MoleScope II is compatible with DermEngine, an intelligent dermatology ecosystem designed for the capturing, documentation and reviewing of skin conditions, including skin cancer.

Please contact us at info@metaoptima.com to request a demo and explore how it will facilitate the improvement of your practice.

You can add and manage your patients, capture and track their skin images, and build your own medical network while providing telemedicine and second opinion services.



Examples of high quality MoleScope II images

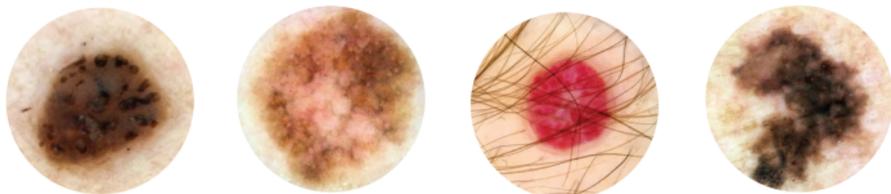

MetaOptima Technology Inc.
Suite 1050, 1185 W Georgia Street,
Vancouver, BC, Canada V6E 4E6
www.metaoptima.com / Tel: +1-778-328-194

EC REP

Emergo Europe
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands



Australia Sponsor
MetaOptima Australia Pty. Ltd.
Level 3, 223 Liverpool Street,
Darlinghurst, NSW, 2010, Australia

UK
CA

UK Responsible Person
Emergo Consulting (UK) Limited,
c/o Cr360 – UL International,
Compass House, Vision Park Histon,
Cambridge CB24 9BZ, United Kingdom