



ALEX® Apps – Dynamic Launch Icons

MetaOptima – DermEngine

Release Date: May 2024

medtechglobal.com

Level 4, 12 Madden Street
Auckland Central
Auckland 1010, New Zealand

PO Box 3329
Shortland Street
Auckland 1140

P 0800 263 3832
P +64 9 358 1123
E support@medtechglobal.com

Contents


Installation Pre-requisite Checklist.....	3
Dynamic Launch Icon Installation Process	4
ALEX® Apps Dynamic Icon Configuration	6
ALEX® Apps Toolbar	7

Installation Pre-requisite Checklist


	Action	Checked
1.	Check the current Medtech Evolution installed version [Help ► About Medtech Evolution]. The practice must be on the following Medtech Evolution version: • Version 6.3 (or higher)	☐
2.	Users are not required to be logged out before running this release. However, users will need to log out and then login for these changes to take effect.	☐
3.	Ensure you are logged onto Windows with Full Administrative rights when performing ANY installation, update or maintenance tasks.	☐
4.	Ensure you are logged onto the Medtech Evolution SQL Server.	☐
5.	Check that this document has been reviewed and understood.	☐

Dynamic Launch Icon Installation Process


1. Ensure that you have completed and checked all the items in the **Installation Prerequisites Checklist** above.
2. Download the **Dynamic Icon Installer-DermEngine.exe** installation file from the location provided.
3. Double-click the **EXE** to start the installation process.
4. If any security warnings appear in your environment, select '**Run**'.
5. The **Dynamic Icon Installer – DermEngine** Installation Wizard will appear.

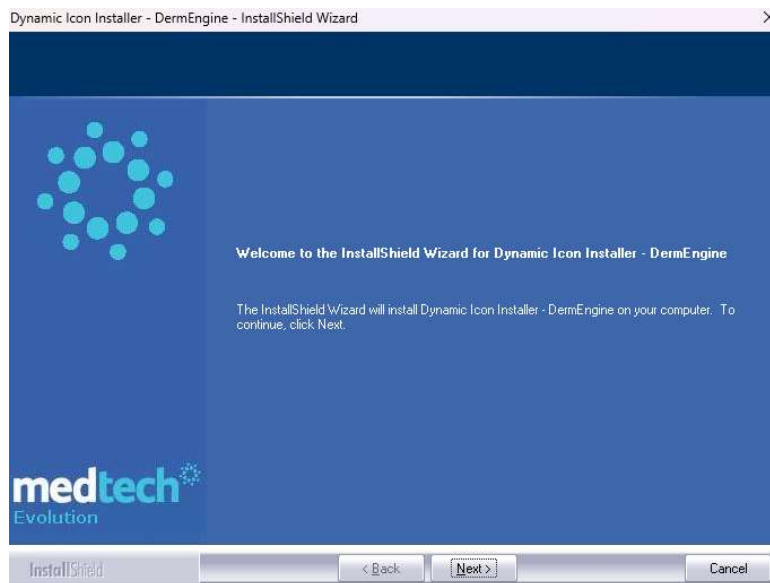

6. Click the '**Next**' button to proceed; the 'Ready to Install' window is displayed.

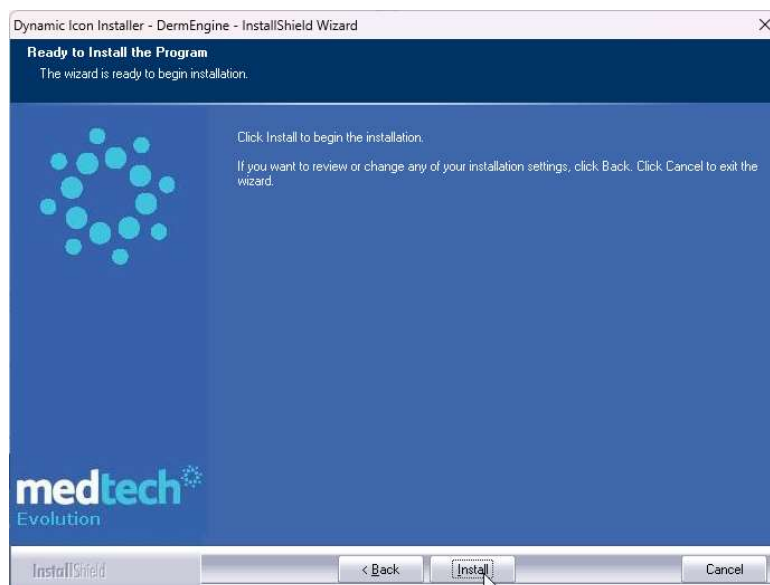


- Click the **'Install'** button to initiate the installation. Once finished the 'Completed' window is displayed.

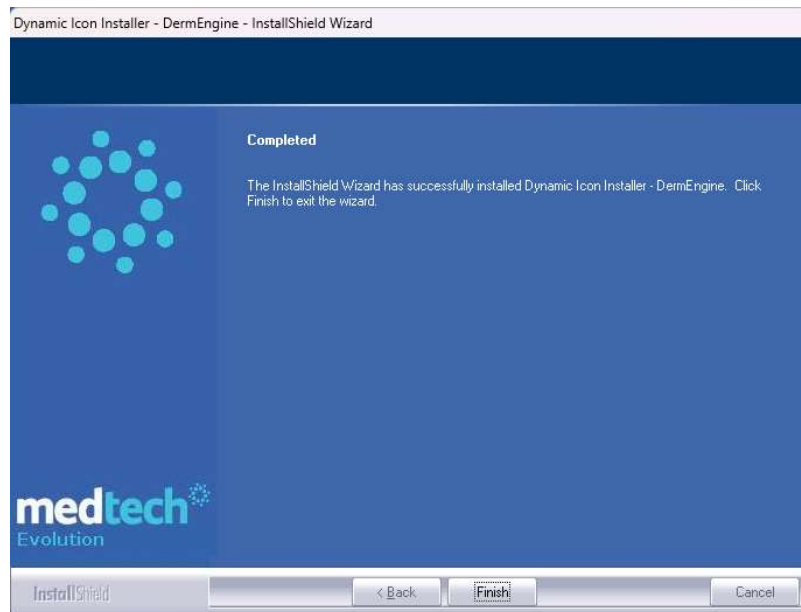

- Click the **'Finish'** button to close the Installation Wizard.

ALEX® Apps Dynamic Icon Configuration

1. Navigate to Evolution icon > Options > Reference Nos > Reference Nos > ALEX® Apps tab.
2. The configuration for the DermEngine icon should be visible in the ALEX® Apps configuration screen.

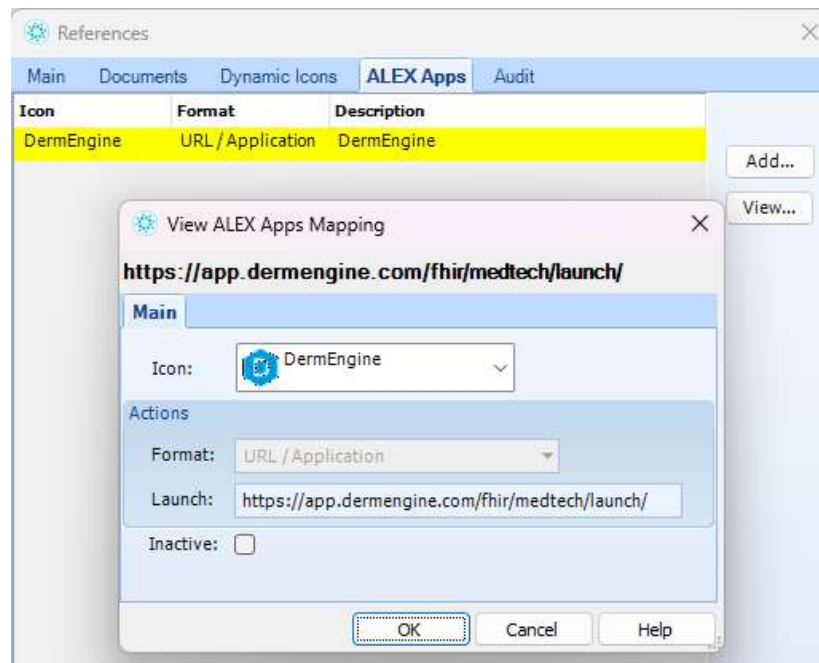


ALEX® Apps Toolbar

1. Ensure that you have access to the ALEX® Apps > Applications option in the Role Based Access Control configuration – this will provide you with access to the ALEX® Apps Toolbar in Evolution.

Role Based Access Control configuration can be found under Evolution icon > Options > Access > Role.

2. The DermEngine icon should be visible in the ALEX® Apps Toolbar for use.

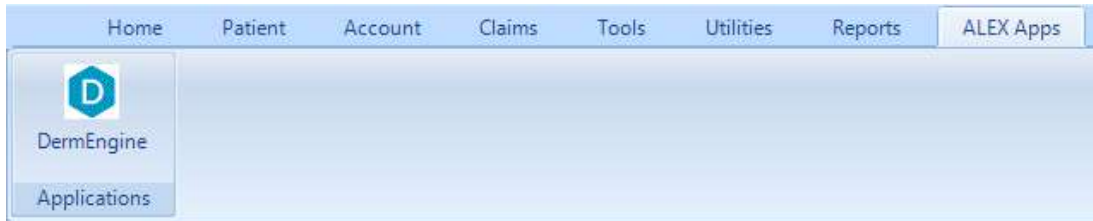

3. Right click on any item in the ALEX® Apps Toolbar and select the Add to Quick Access Toolbar option to add the selected icon onto your user specific Quick Access Toolbar.

